

#### Specifications

Power Supply	100~240VAC, 50/60Hz, 200VA
Flow Rate	0.1~10ml/s
Pressure Limit	325PSI
Auto Filling Rate	8ml/s
Air Purge Rate	1.0~9.9ml/s
Piston Rod Speed	9.9ml/s
Phases	6
Injection Programs	2000
Injection Records	2000
Volume	A: 1~200ml
	B: 1~200ml
Multiple Language User Interface	

#### High Quality Disposables



200ml Syringe for Contrast
200ml Syringe for Saline

## antmed

# Shenzhen Antmed Co., Ltd. 18 Jinhui Ave., Pingshan New District, Shenzhen 518122, China

Tel: +86 755 8606 0992 Fax: +86 755 8606 0991 E-mail: info@antmed.com http://www.antmed.com





### CT Dual Injector

Leading Contrast Media Injectors Forward



## CT Dual Injector

· Able to observe the right coronary arteries,

right atrium and right ventricle in the same image

Introducing an all-new ImaStar CT dual injector, offering an extended list of product features, reliable performance, and a safe patient experience.



#### Color Touch-Screen

- 12.1 and 7" color touch-screens for easy operation
- Clear and concise information display to ensure user operational confidence

### Injection Record Storage

· Stores up to 2000 injection history records



#### Auto Detect


· Syringe detection feature automates air purging and piston rod operation

### Connect to Other Imaging Equipment

- Connect with other imaging equipment for enhanced imaging, use contrast more efficiently and minimize radiation
- Option


### Real-time Pressure Curve

 Display provide real-time pressure monitoring curve

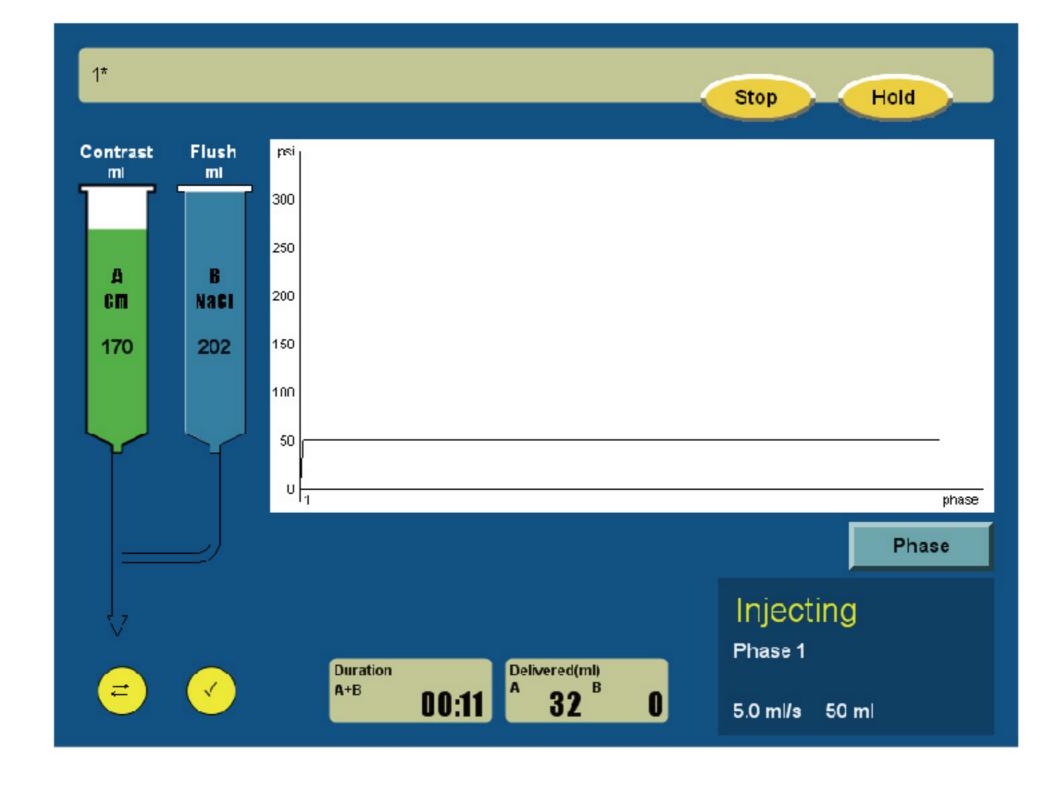

# A Variety of Configuration Options

- Stores up to 2000 customized injection programs
- Program storage capacity can be increased based on user needs

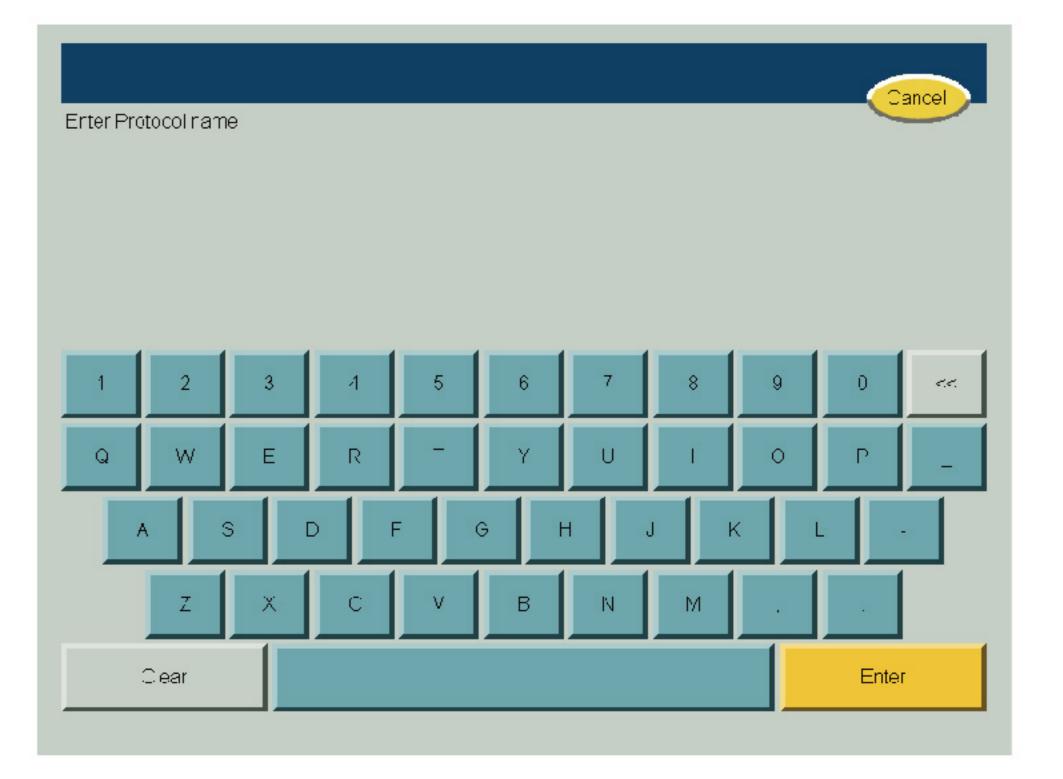

### Safety

- · To ensure the safety of the injection, there is a double confirmation feature with separate prompts to confirm air purging and injection ready status before any injection
- · Automatically slows down the injection when approaching to pressure limit
- · Injection stops when pressure limit is reached. Alerts with both light and sound
- · Waterproof injector head designed to minimize damage from contrast and saline leakage
- · Ability to pause and resume while injection is underway